

# Data from a proteomic analysis highlight different osmoadaptations in two strain of Propionibacterium freudenreichii

Floriane Gaucher, Sylvie Bonnassie, Houem Rabah, Pauline Leverrier, Sandrine Pottier, Julien Jardin, Valérie Briard-Bion, Pierre Marchand, Romain Jeantet, Philippe Blanc, et al.

# ▶ To cite this version:

Floriane Gaucher, Sylvie Bonnassie, Houem Rabah, Pauline Leverrier, Sandrine Pottier, et al.. Data from a proteomic analysis highlight different osmoadaptations in two strain of Propionibacterium freudenreichii. Data in Brief, 2020, 28, pp.104932. 10.1016/j.dib.2019.104932 . hal-02404929

# HAL Id: hal-02404929

# https://institut-agro-rennes-angers.hal.science/hal-02404929

Submitted on 11 Dec 2019

**HAL** is a multi-disciplinary open access archive for the deposit and dissemination of scientific research documents, whether they are published or not. The documents may come from teaching and research institutions in France or abroad, or from public or private research centers. L'archive ouverte pluridisciplinaire **HAL**, est destinée au dépôt et à la diffusion de documents scientifiques de niveau recherche, publiés ou non, émanant des établissements d'enseignement et de recherche français ou étrangers, des laboratoires publics ou privés.



Distributed under a Creative Commons Attribution - NonCommercial - ShareAlike 4.0 International License

Data in brief xxx (xxxx) xxx



Contents lists available at ScienceDirect

# Data in brief

journal homepage: www.elsevier.com/locate/dib

## Data Article

# Data from a proteomic analysis highlight different osmoadaptations in two strain of Propionibacterium freudenreichii

- **6** Floriane Gaucher <sup>a, b</sup>, Sylvie Bonnassie <sup>a, c</sup>, Houem Rabah <sup>a, d</sup>, Pauline Leverrier <sup>e</sup>, Sandrine Pottier <sup>f</sup>, Julien Jardin <sup>a</sup>, Valérie Briard-Bion <sup>a</sup>, Pierre Marchand <sup>b</sup>, Romain Jeantet <sup>a</sup>, Philippe Blanc <sup>b</sup>, Gwénaël Jan <sup>a, \*</sup>
  - <sup>a</sup> UMR STLO, Agrocampus Ouest, INRA, F-35042, Rennes, France
  - <sup>b</sup> Bioprox, 6 rue Barbès, 92532, Levallois-Perret, France
  - <sup>c</sup> Université de Rennes I, Univ. Rennes, Rennes, France
  - <sup>d</sup> Bba, Pôle Agronomique Ouest, Régions Bretagne et Pays de la Loire, F-35042 Rennes, France
  - <sup>e</sup> de Duve Institute, Université catholique de Louvain, Avenue Hippocrate 75, Brussels, 1200, Belgium
  - <sup>f</sup> Univ. Rennes, CNRS, ISCR UMR 6226, PRISM, BIOSIT UMS 3480, F-35000, Rennes, France

#### ARTICLE INFO

Article hist	ory:
Received 3	30 August 2019
Received in	n revised form 26 November 2019
Accepted 2	27 November 2019
Available o	online xxx

Keywords:
Propionibacterium
Osmotic
Stress
Probiotic
Osmoadaptation
Proteomic
Cheese
Osmoprotectant
-

#### ABSTRACT

The article presents a proteomic data set generated by a comparative analysis of the proteomes of *Propionibacterium freudenreichii*, comparing the CIRM-BIA 129 and CIRM-BIA 1025 strains. The two strains were cultivated until the beginning of the stationary phase in a chemical defined medium (MMO), and in this medium in the presence of NaCl, with or without glycine betaine. Whole-cell proteins were extracted, trypsinolyzed and analyzed by nano LC-MS/MS, prior to identification and classification by function using the X!Tandem pipeline software and the proteomic data from NCBI.nlm.nigh.gov. Quantification of proteins was then carried out in order to detect change in their expression depending on the culture medium. This article is related to the research article entitled "Benefits and drawbacks of osmotic adjustment in *Propionibacterium freudenreichii*". The comparative proteomic analysis

\* Corresponding author.

E-mail address: gwenael.jan@inra.fr (G. Jan).

#### https://doi.org/10.1016/j.dib.2019.104932

2352-3409/© 2019 Published by Elsevier Inc. This is an open access article under the CC BY-NC-ND license (http:// creativecommons.org/licenses/by-nc-nd/4.0/).

Please cite this article as: F. Gaucher et al., Data from a proteomic analysis highlight different osmoadaptations in two strain of *Propionibacterium freudenreichii*, Data in brief, https://doi.org/10.1016/ j.dib.2019.104932

42

43

F. Gaucher et al. / Data in brief xxx (xxxx) xxx

of the two strains reveal strain-dependent and mediumdependent stress proteomes in the probiotic *P. freudenreichii*. © 2019 Published by Elsevier Inc. This is an open access article under the CC BY-NC-ND license (http://creativecommons.org/ licenses/by-nc-nd/4.0/).

Specifications Table

Subject area	Microbiology, Food technology
More specific subject area	Probiotics, dairy starters
Type of data	Table and figure
How data was acquired	Mass spectrometry
Data format	Raw, Analyzed and Filtrated
Experimental factors	Propionibacterium freudenreichii CIRM-BIA 129 and CIRM-BIA 1025 were grown in chemically defined medium, and in presence of salt with or without glycine betaine. Whole cell proteins were extracted and analyzed by nanoLC-MS/MS in a gel-free strategy.
Experimental features	All the samples were analyzed by nano-LC coupled to MS/MS
Data source location	Both strains of <i>P. freudenreichii</i> , isolated from French cheeses, are the property of CNIEL, France. Data were acquired at INRA STLO, France, under the agreement of an MTA.
Data accessibility	Data have been deposit into an INRA data base: https://doi.org/10.15454/VVNQHR
Related research article	Gaucher F., Bonnassie S, Rabah H, Leverrier P, Pottier, S, Jardin J, Briard-Brion V, Marchand P, Jeantet R, Blanc, P, Jan G, Benefits and drawbacks of osmotic adjustment in <i>Propionibacterium freudenreichii</i> , Journal of proteomics (2019). In press

#### Value of the Data

• This dataset is an important step in understanding the molecular mechanisms responsible for osmoadaptation in P. freudenreichii.

• A chemically defined medium was developed to clearly distinguish the effect of 1) osmotic constraint imposed by salt and 2) osmoadaptation in the presence of the universal osmoprotectants glycine betaine in a probiotic bacterium.

• Surprisingly, molecular mechanisms of osmoadaptation are shown here to vary among strains within the same species.

#### 1. Data

This article presents a dataset generated during a proteomic comparative analysis of the cellular proteome of two phenotypical different *Propionibacterium freudenreichii* strains (CIRM-BIA 129 and CIRM-BIA 1025). The list of all proteins identified is provided in Supplementary Table 1, which presents the protein relative abundances. Fig. 1 shows a heatmap corresponding to abundances, for the two strains and the three culture media (chemically defined medium, the same with salt, the same with salt and glycine betaine). This report presents data from the research article entitled "Benefits and drawbacks of osmotic adjustment in *Propionibacterium freudenreichii*". The cellular proteomes reveal diversity in the osmoadaptation of *P. freudenreichii*, among strains.

The proteomic data acquired during this study are presented in Supplementary Table 1 which summarizes all the mean values and all protein with a ratio above 1.5 or below 0.66. For each biological and each technical replicate, a raw data file is presented (https://doi.org/10.15454/VVNQHR).


Fig. 1 summarizes all the observed inductions and repressions. It is composed of a heatmap corresponding to abundances, for the two strains and the three culture media (chemically defined medium, the same with salt, the same with salt and glycine betaine).

Please cite this article as: F. Gaucher et al., Data from a proteomic analysis highlight different osmoadaptations in two strain of *Propionibacterium freudenreichii*, Data in brief, https://doi.org/10.1016/ j.dib.2019.104932

## ARTICLE IN PRESS DIB104932\_proof = 10 December 2019 = 3/5


#### F. Gaucher et al. / Data in brief xxx (xxxx) xxx

3

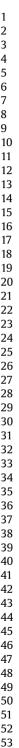

Q4

54

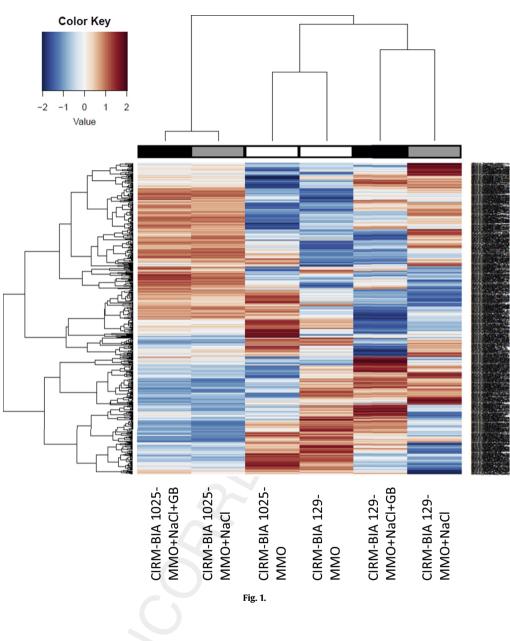

### 2. Experimental design, materials, and methods

### 2.1. Strain and pre-culture

Propionibacterium freudenreichii CIRM-BIA 129 (equivalent ITG P20) and CIRM-BIA 1025 (equivalent ITG P1) were provided, stored and maintained by the CIRM-BIA Biological Resource Center (Centre International de Ressources Microbiennes-Bactéries d'Intérêt Alimentaire, INRA, Rennes, France). The strains are routinely cultivated in yeast extract lactate (YEL) medium [1] in this study, they are also

Please cite this article as: F. Gaucher et al., Data from a proteomic analysis highlight different osmoadaptations in two strain of *Propionibacterium freudenreichii*, Data in brief, https://doi.org/10.1016/ j.dib.2019.104932 

#### F. Gaucher et al. / Data in brief xxx (xxxx) xxx

cultivated in MMO (Medium Minus Osmoprotectants). MMO is a synthetic medium which was derivated from CdM (Chemiclal defined Medium) describe previously [2], MMO medium has the same composition as CdM but glutamate, glutamine and proline, which are considered as potent sources of osmoprotectants, were removed. *P. freudenreichii* was grown at 30 °C without agitation under microphilic condition. Composition of MMO is detailed in Supplementary Table 1.

#### 2.2. Whole cell protein extraction and trypsin digestion

The label free proteomics has have been conducted as Huang et al., 2018 [3]. At the beginning of stationary phase, *P. freudenreichii* cells were harvested by centrifugation and washed twice with PBS buffer (NaCl 8g,L<sup>-1</sup>, KCl 2g,L<sup>-1</sup> KH<sub>2</sub>PO<sub>4</sub> 2g,L<sup>-1</sup>, Na<sub>2</sub>HPO<sub>4</sub> 12H<sub>2</sub>O 35,8g,L<sup>-1</sup>). Cell pellets were then resuspended in lysis solution (lysis solution: 0.5 mL pH 7.5, 157g Tris-HCl, 0.03 g SDS, 0.3 g DTT and 9.5 mL H<sub>2</sub>O), with 1 mM phenylmethylsulfonyl fluoride (PMSF, Sigma-Aldrich, USA). The solution was frozen for 1 hour, then sonicated (2 min 30 HZ), and cells were broken using zirconium beads in the homogenizer (Homogénéisateur Precellys Evolution - Bertin Instruments, France). The resulting SDS extracts were recovered by centrifugation (21,000 ×g; 20 °C; 20 min) and then cleaned and quantified using the two-dimensional (2-D) clean-up kit (GE Healthcare Bio-Sciences AB, Uppsala, Sweden) and the 2-D quant kit (GE Healthcare Bio-Sciences AB, Uppsala, Sweden) and the 2-D quant kit (GE Healthcare Bio-Sciences AB, Uppsala, Sweden), respectively. Tryptic digestion was performed on 100 µg of whole-cell proteins from each sample during 15 hours at 37 °C using Sequencing Grade Modified Trypsin (Promega, Madison, USA) according to the manufacturer's instructions and as described previously [3]. Spectrophotometric-grade trifluoroacetic acid (TFA)(Sigma-Aldrich, USA) was added in order to stop the digestion.

#### 2.3. Nano LC-MS/MS

Experiments were performed as previously described [3]. Experiments were performed using a nano RSLC Dionex U3000 system fitted to a Q-Exactive mass spectrometer (Thermo Scientific, San Jose, USA) equipped with a nano-electrospray ion source. A preliminary sample concentration step was performed on a C18 pepMap100 reverse phase column (C18 column, 300-µm inner diameter (i.d.), 5 mm length, 5 µm particle size, 100 Å pore size; Dionex, Amsterdam, The Netherlands). Peptides separation was performed on a reversed-phase column (PepMap 100 C18, 75 μm i.d., 250 mm length, 3 μm particle size, 100 Å pore size; Dionex, Amsterdam, The Netherlands) with a column temperature of 35 °C, using solvent A (2% (v/v) acetonitrile (Honeywell, USA), 0.08% (v/v) formic acid and 0.01% (v/v) TFA in deionized water) and solvent B (95% (v/v) acetonitrile, 0.08% (v/v) formic acid and 0.01% (v/v) TFA in deionized water). Peptides were separated using a gradient of 5–35% solvent B over 80 min followed by 35-85% solvent B over 5 min at a flow rate of 0.3  $\mu$ L/min. Eluted peptides were directly electro sprayed into the mass spectrometer operating in positive ion mode with a voltage of 2kV. Spectra were recorded in full MS mode and selected in a mass range 250-2000 m/z for MS spectra with a resolution of 70,000 at m/z 200. For each scan, the ten most intense ions were selected for fragmentation. MS/MS spectra were recorded with a resolution of 17,500 at m/z 200 and the parent ion was subsequently excluded from MS/MS fragmentation for 20 s. The instrument was externally calibrated according to the supplier's instructions.

### 2.4. Proteins identification

Proteins identification was performed as previously described [3]. Peptides were identified from the MS/MS spectra using the X!Tandem pipeline software (Langella et al., 2017). The search was performed against a database composed of proteomes of strains P. freudenreichii CIRM-BIA 129 and CIRM-BIA 1025 (downloaded from NCBI.nlm.nih.gov on the August 23, 2018). Database search parameters were specified as follow: trypsin cleavage was used and the peptide mass tolerance was set to 10 ppm for MS and 0.05 Da for MS/MS. Oxidation of methionine and phosphorylation of threonine, serine and tryptophan were selected as a variable modification. For each peptide identified, a minimum score corresponding to an e-value below 0.05 was considered as a prerequisite for peptide validation.

Please cite this article as: F. Gaucher et al., Data from a proteomic analysis highlight different osmoadaptations in two strain of *Propionibacterium freudenreichii*, Data in brief, https://doi.org/10.1016/ j.dib.2019.104932 F. Gaucher et al. / Data in brief xxx (xxxx) xxx

#### 2.5. Protein quantification

1

2 3

4

5

6

7

8 9

10

11

12

13

14

15

16 17

18

19

20 21

22 23

24 25

26

27 28

29 30

31 32

33

34 35

36 37

38 39

40

41

42

43

44

45

46

47

48

49

50 51 52 02

Protein quantification was performed as previously described [3]. Every peptide identified by tandem mass spectrometry was quantified using the free MassChroQ software (Valot et al., 2011) before data treatment and statistical analysis within the R software (R 3.2.2, Project for statistical computing). A specific R package called 'MassChrogR' was used to automatically filter dubious peptides for which standard deviation of retention time was superior to 30 s and to regroup peptide quantification data into proteins. Two different and complementary methods of analysis were used, based on peak counting or XIC (eXtracted Ion Current). For peak counting analysis, variance analysis was performed to detect significant abundance variation between culture conditions. Proteins with an adjusted p-value <0.05 were considered significantly different. For XIC based quantification, normalization was performed to take into account possible global quantitative variations between LC-MS runs. Peptides shared between different proteins were automatically excluded from the data set as well as peptides present in less than 95% of samples. Missing data were then imputed from a linear regression based on other peptide intensities for the same protein [4]. Protein abundance was computed as the sum of peptide areas requiring a minimum of 2 peptides per protein. For every missing protein abundance, missing value were replaced by the minimum abundance obtained for this protein in the whole experiment. Analysis of variance was used to determine proteins with significantly different abundance between our culture conditions.

### Uncited reference

## **Q5** [5].

#### Acknowledgments

F. Gaucher is the recipient of a joint doctoral fellowship from Bioprox and from ANRT (Association Nationale de la Recherche et de la Technologie).

### **Conflict of interest**

The authors declare that they have no known competing financial interests or personal relationships that could have appeared to influence the work reported in this paper.

### Appendix A. Supplementary data

Supplementary data to this article can be found online at https://doi.org/10.1016/j.dib.2019.104932.

#### References

- A.C. Malik, G.W. Reinbold, E.R. Vedamuthu, An evaluation of the taxonomy of Propionibacterium, Can. J. Microbiol. 14 (1968) 1185–1191.
- [2] V. Gagnaire, J. Jardin, H. Rabah, V. Briard-Bion, G. Jan, Emmental cheese environment enhances Propionibacterium freudenreichii stress tolerance, PLoS One 10 (2015) e0135780, https://doi.org/10.1371/journal.pone.0135780.
- [3] S. Huang, F. Gaucher, C. Cauty, J. Jardin, Y. Le Loir, R. Jeantet, X.D. Chen, G. Jan, Growth in hyper-concentrated sweet whey triggers multi stress tolerance and spray drying survival in Lactobacillus casei BL23: from the molecular basis to new perspectives for sustainable probiotic production, Front. Microbiol. (2018).
- [4] M. Blein-Nicolas, W. Albertin, T. da Silva, B. Valot, T. Balliau, I. Masneuf-Pomarède, M. Bely, P. Marullo, D. Sicard, C. Dillmann, D. de Vienne, M. Zivy, A systems approach to elucidate heterosis of protein abundances in yeast, Mol. Cell. Proteom. 14 (2015) 2056–2071, https://doi.org/10.1074/mcp.M115.048058.
- [5] F. Gaucher, S. Bonnassie, H. Rabah, P. Leverrier, S. Pottier, J. Jardin, V. Braird-Brion, P. Marchand, R. Jeantet, P. Blanc, G. Jan, Benefits and drawbacks of osmotic adjustment in *Propionibacterium freudenreichii*, Journal of proteomics (2019). In press.